Tru-Incise[®] Valvulotome Cutting Head Kit

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Tru-Incise® Valvulotome Cutting Head Kit (Model # TIVK2030)

Instructions For Use - English

Rx only **STERILE** **EO** 


Introduction

Preparation of the saphenous vein for an in situ arterial bypass graft requires the disruption of the valve sets with minimal damage to the intima of the vein. The Tru-Incise Valvulotome is designed with this goal in mind. Features of the Tru-Incise Valvulotome include the irrigation port at the end of the cutter, interchangeable cutting heads designed to optimize the cutting surface relative to a tapering vein, and the pronged blades that engage the valve sets for cutting, thus maximizing cutting action.

Accessories Include:

- * 1 Valvulotome Catheter with vein sleeve (94cm)
- * 1 Blunt introducing head
- * 1 2 mm cutting head
- * 1 3 mm cutting head
- * 1 4 mm cutting head
- * 1 2 mm tributary (head introducer) sleeve
- * 1 3 mm tributary (head introducer) sleeve
- * 1 4 mm tributary (head introducer) sleeve

Valve disruption can be performed prior to or following anastomosis.

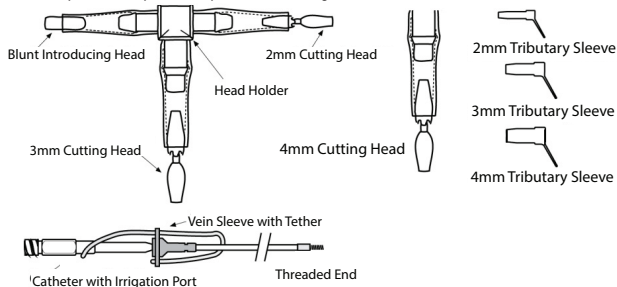The Tru-Incise Valvulotome Catheter Kit (REF # TIVK2030) provides the flexibility of interchangeable cutting heads to allow appropriate sizing for varying vessel diameters. These instructions for use pertain to use of the catheter kit which includes individual cutting heads.

Sterility

Sterile if package is unopened and undamaged. Non-Pyrogenic. For single use only. The reuse of this single-use device can lead to patient infection and/or device malfunction.

Indication

In an in situ bypass procedure, the valves must be rendered incompetent to allow retrograde flow of arterial blood in the vein. The Tru-Incise Valvulotome is specifically indicated for use during operations requiring rendering valve leaflets incompetent while minimizing intimal damage. This is primarily intended for use during in situ bypass procedures and may be utilized during distal infrainguinal bypass when a non-anatomic position is required (ex: profunda to anterior tibial artery), composite vein infrainguinal bypass, or even an aorto-renal bypass.

Contraindications

This device is not designed to be used for endarterectomy, vein stripping, embolectomy, vessel dilation, or thrombolysis.

To Open the Package

Grasping the outer pouch, peel the pouch open exposing the tray. Using sterile technique, carefully remove the tray while placing it in a sterile field. Open the tray by peeling back the tray lid at the tray tab, exposing the contents. Carefully remove the catheter and the remaining contents as needed.


Possible Complications with the Use of the Valvulotome

Complications which may occur with the use of any valvulotome include, but are not limited to: intimal disruption; vessel wall perforation; missed/incomplete valve disruption; phlebitis; thrombosis; stricture; A-V fistulae; entrapment of cutting heads; infection.


SUGGESTED PROCEDURE - A: VALVE DISRUPTION PRIOR TO ANASTOMOSIS (Procedure)

This procedure relies on the use of delivered irrigant to ensure closure of venous valve sets prior to disruption. It should be noted that these instructions are not limited to femoral-to-below-knee bypasses. The techniques of manipulation of the valvulotome for rendering valves incompetent are equally applicable to other bypasses requiring valve disruption.

1. Expose and prepare the saphenous vein using standard technique. If the Tru-Incise Valvulotome is not going to be used to disrupt the first valve set, cut the first valve set under direct vision using standard technique.
2. Using the head holder, attach the blunt introducing head to the end of the catheter. Fully tighten the head prior to use (approximately eight full rotations).
3. Introduce the catheter with the blunt introducing head through the below knee end of the prepared vein. If desired, the vein sleeve can be affixed to the cut end of the vein using the attached tether to assist in easy passage of the catheter through the open end of the vein. Irrigate with heparinized saline (or other suitable solution) through the catheter irrigation port to dilate the vein prior to advancing the device into the vein.

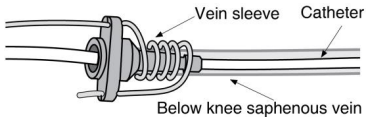
4. Gently clamp the end of the vein severed at the fossa ovalis from the deep vein to allow for vein dilation. With a 20 cc syringe, use gentle irrigation to minimize hydrostatic pressure injuries to the vessel wall. With the vein distended, gently advance the catheter while allowing it to float up the main channel of the saphenous vein.

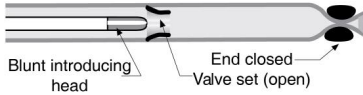
5. When the catheter reaches the clamped sapheno-femoral vein end, release the clamp, advance the blunt introducing head beyond the end of the cut vein, and replace the blunt introducing head with an appropriately sized valve cutter head. To minimize the risk of cutting through a glove and to ensure proper torque, the head holder should be used to attach the cutting head. Fully tighten the cutting head prior to use (approximately eight full rotations). DO NOT ROTATE THE CATHETER DURING USE. Rotation could cause unscrewing of the cutting head and/or damage to the vessel wall.

NOTE: The surgeon shall determine the appropriately sized cutting head which will best cut the valve sets at a given position in the vessel with the least damage to the vessel wall.

6. Position the valve cutter at the saphenofemoral end of the vein while injecting fluid through the catheter. This will distend the lumen and close the next valve set which will now be appropriately positioned for cutting.

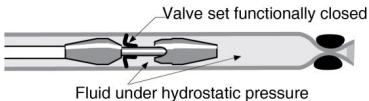7. Pull the valve cutter back to engage and cut the first or next valve set. Ensure that the section of vein leading up to that valve set is fully distended to close the valves. Use slow and consistent traction to disrupt valve sets. While maintaining the hydrostatic pressure, pull the valve cutter to engage and cut each sequential valve set until all valves have been rendered incompetent within the appropriate section of vein relative to the chosen cutting head size.
8. Change cutting heads as needed during the procedure. Cutting heads may be exchanged through the saphenofemoral end of the vein or through a major tributary using the appropriate tributary sleeve. Appropriately match the size of the cutting head to the vessel's tapering lumen. The clearance between the vein and cutter head should be approximately 1/2 mm.

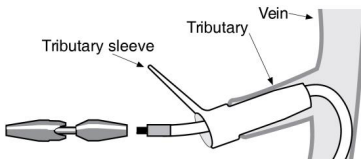9. Withdraw the valve cutter until it engages the vein sleeve and remove the engaged cutter/sleeve assembly.


10. Ligate the venous tributaries.
11. When all valve sets are rendered incompetent, confirm free flow of irrigant through the saphenous vein. Complete the below knee anastomosis to the appropriate arterial segment. Confirm valve disruption and tributary ligation by doppler ultrasound, fluoroangiography, angiography and/or duplex doppler.
12. Interrupt any significant remaining A-V fistulae to ensure adequate bypass flow.

SUGGESTED PROCEDURE - B: VALVE DISRUPTION FOLLOWING ANASTOMOSIS


This procedure relies on systemic blood pressure to ensure venous valve set closure prior to disruption.

It should be noted that these instructions are not limited to femoral-to-below-knee bypasses. The techniques of manipulation of the valvulotome for rendering valves incompetent are equally applicable to other bypasses requiring valve disruption.

1. Expose and prepare the saphenous vein using standard technique. If the Tru-Incise Valvulotome is not going to be used to disrupt the first valve set, cut the first valve set under direct vision using standard technique.
2. Using the head holder, attach the blunt introducing head to the end of the catheter. Fully tighten the head prior to use (approximately eight full rotations).
3. Introduce the catheter with the blunt introducing head through the below knee end of the prepared vein. If desired, the vein sleeve can be fixed to the cut end of the vein using the attached tether to assist in easy passage of the catheter through the open end of the vein. Irrigate with heparinized saline (or other suitable solution) through the catheter irrigation port to dilate the vein prior to advancing the device through the vein.

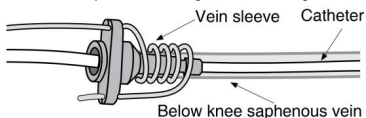
4. Gently clamp the end of the vein severed at the fossa ovalis from the deep vein to allow for vein dilation. With a 20 cc syringe, use gentle irrigation to minimize hydrostatic pressure injuries to the vessel wall. With the vein distended, gently advance the catheter while allowing it to float up the main channel of the saphenous vein.

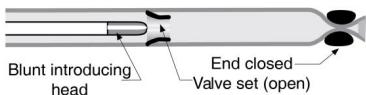
5. When the catheter reaches the clamped sapheno-femoral vein end, release the clamp, advance the blunt introducing head beyond the end of the cut vein, and replace the blunt introducing head with an appropriately sized valve cutter head. To minimize the risk of cutting through a glove and to ensure proper torque, the head holder should be used to attach the cutting head. Fully tighten the cutting head prior to use (approximately eight full rotations). DO NOT ROTATE THE CATHETER DURING USE. Rotation could cause unscrewing of the cutting head and/or damage to the vessel wall.

NOTE: The surgeon shall determine the appropriately sized cutting head which will best cut the valve sets at a given position in the vessel with the least damage to the vessel wall.

6. Perform an anastomosis of the open saphenofemoral vein end to the common femoral or other appropriate artery. Blood flow through the anastomosis will be retained at the first valve encountered.
7. Withdraw the catheter to render the first or next valve set incompetent. Use slow and consistent traction to disrupt valve sets. Valve disruption will be confirmed as blood flow will now be retained at the site of the next restraining valve.

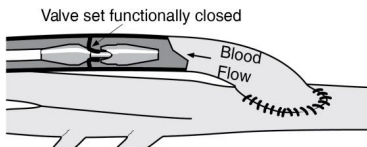8. Change cutting heads as needed during the procedure. Cutting heads may be exchanged through a major tributary using the appropriate tributary sleeve or through a partially completed saphenofemoral anastomosis. Appropriately match the size of the cutting head to the vessel's tapering lumen. The clearance between the vein and cutter head should be approximately 1/2 mm.

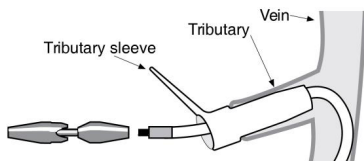9. Disrupt each valve set in a like manner.
10. Withdraw the valve cutter until it engages the vein sleeve and remove the engaged cutter/sleeve assembly.
11. Ligate the venous tributaries.
12. When all valve sets are rendered incompetent, confirm free flow of irrigant through the saphenous vein. Complete the below knee anastomosis to the appropriate arterial segment. Confirm valve disruption and tributary ligation by doppler ultrasound, fluoroangiography, angiography and/or duplex doppler.
13. Interrupt any significant remaining A-V fistulae to ensure adequate bypass flow.

WARNINGS/PRECAUTIONS

1. Using an inappropriately sized cutting head may result in damage to the vein wall.
2. Veins are known to change diameter significantly in some individuals. Determine and consider the patient's anatomy when selecting valve cutter size.
3. During valve disruption, gently withdraw the catheter.
4. Do not over-distend the vessel with irrigation.
5. Fully tighten the cutting head prior to use. DO NOT ROTATE THE CATHETER DURING USE. Rotation could cause unscrewing of the cutting head and/or damage to the vessel wall.
6. Use the head holder to minimize the risk of cutting through a glove.
7. When using more than one cutter head to lyse valves, initiate the procedure with the largest head. This will aid in minimizing inadequate incising of valve leaflets.
8. When removing heads from the head holder, always grasp the catheter in one hand and the CENTER of the rigid portion of the head holder in the other and pull straight out.

Resterilization/Re-use

This device is single-use only. Do not reuse, reprocess, or re-sterilize. The cleanliness and sterility of the re-processed device cannot be assured. Reuse of the device may lead to cross contamination, infection, or patient death. The performance characteristics of the device may be compromised due to reprocessing or re-sterilization since the device was only designed and tested for single use. The shelf life of the device is based on single use only.

Notices: Limited Product Warranty; Limitation Of Remedies

LeMaitre Vascular, Inc. warrants that reasonable care has been used in the manufacture of this device. Except as explicitly provided herein, LEMAITRE VASCULAR (AS USED IN THIS SECTION, SUCH TERM INCLUDES LEMAITRE VASCULAR, INC., ITS AFFILIATES, AND THEIR RESPECTIVE EMPLOYEES, OFFICERS, DIRECTORS, MANAGERS, AND AGENTS) MAKES NO EXPRESS OR IMPLIED WARRANTIES WITH RESPECT TO THIS DEVICE, WHETHER ARISING BY OPERATION OF LAW OR OTHERWISE (INCLUDING, WITHOUT LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE) AND HEREBY DISCLAIMS THE SAME. LeMaitre Vascular makes no representation regarding the suitability for any particular treatment in which this device is used, which determination is the sole responsibility of the purchaser. This limited warranty does not apply to the extent of any abuse or misuse of, or failure to properly store, this device by the purchaser or any third party. The sole remedy for a breach of this limited warranty shall be replacement of, or refund of the purchase price for, this device (at LeMaitre Vascular's sole option) following the purchaser's return of the device to LeMaitre Vascular. This warranty shall terminate on the expiration date for this device.

IN NO EVENT SHALL LEMAITRE VASCULAR BE LIABLE FOR ANY DIRECT, INDIRECT, CONSEQUENTIAL, SPECIAL, PUNITIVE, OR EXEMPLARY DAMAGES. IN NO EVENT WILL THE AGGREGATE LIABILITY OF LEMAITRE VASCULAR WITH RESPECT TO THIS DEVICE, HOWEVER ARISING, UNDER ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, TORT, STRICT LIABILITY, OR OTHERWISE, EXCEED ONE THOUSAND DOLLARS (US\$1,000), REGARDLESS OF WHETHER LEMAITRE VASCULAR HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH LOSS, AND NOTWITHSTANDING THE FAILURE OF THE ESSENTIAL PURPOSE OF ANY REMEDY. THESE LIMITATIONS APPLY TO ANY THIRD-PARTY CLAIMS.

A revision or issue date for these instructions is included on the back page of these Instructions for Use for the user's information. If twenty-four (24) months has elapsed between this date and product use, the user should contact LeMaitre Vascular to see if additional product information is available.

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LeMaitre and Tru-Incise are registered trademarks of LeMaitre Vascular, Inc.
©2016 LeMaitre Vascular, Inc. All rights reserved. Printed in the USA.
LeMaitre Vascular GmbH