**AnastoClip GC<sup>®</sup> Closure System**

English — Instructions for Use

**AnastoClip GC<sup>®</sup> Verschlussystem**

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**STERILE EO Rx only**

## IMPORTANT!

This booklet is designed to assist in using the AnastoClip GC® Closure System with titanium clips. It is not a reference to surgical stapling techniques.

BEFORE USING PRODUCT, READ THE FOLLOWING INFORMATION THOROUGHLY.

## Indications


The AnastoClip GC is intended for use in the closure of arteriotomies and veinotomies, the attachment of synthetic vascular prostheses, and the creation of everting anastomoses in blood vessels and other small tubular structures. The Applier is also intended for approximation of dural tissue (dural approximation only available in the European Union).

## Effects

The AnastoClip GC applier is available in three (3) clip sizes: medium-1.1 mm, large-1.7 mm and extra large-2.5 mm. They are available in two (2) lengths: 3 in. distal shaft and 6 in. distal shaft. The AnastoClip GC applier consists of a rotating shaft and an integral cartridge containing titanium clips. The AnastoClip implantable clip is made of titanium and it is MRI-Conditional (up to 3 Tesla). As the levers of the applier are squeezed together, the clip is closed around the everted tissue edges. As the levers are released, a new clip is automatically loaded into the clip applier jaws. It is recommended to use the Tissue Everting Forceps to aid in the everting of the tissue edges and the AnastoClip® Remover for the removal of any AnastoClip GC clips (if necessary).

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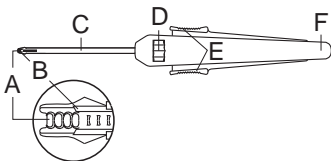## How Supplied:

The applier is supplied sterile. The sterility of the device is assured as long as the packaging is not opened or damaged.

The entire device is considered to be non-pyrogenic for vascular applications.

The implantable clips, the stainless steel shaft and jaws of the applier are considered to be non-pyrogenic for dura applications. The plastic part of the applier (handle, levers and rotation knob) should not be in contact with cerebral spinal fluid.

## Instructions For Use


**NOTE:** It is recommended to wear loupes. A 2.5X magnification is suggested.

1. Preparation of tissues is recommended as follows:

- G) ARTERIOTOMY OR VENOTOMY: One optional stay suture at mid-incision.
- H) END-TO-END: Horizontal mattress sutures at 3 and 9 o'clock.
- I) END-TO-SIDE: Horizontal mattress sutures heel and toe: stay sutures at 3 and 9 o'clock.
- J) SIDE-TO-SIDE: Horizontal mattress sutures at 12 and 6 o'clock: stay sutures at 3 and 9 o'clock.

**NOTE:** Additional sutures may be placed depending on the length of the closure, in order to facilitate a symmetrical eversion.


- 2. Symmetrically evert all tissue layers for secure nonpenetrating clip placement. Evert the tissue edges of the vessel with either design of tissue everting forceps. Ensure that all tissue edges are symmetrically everted prior to applying the clip. Failure to symmetrically evert the tissue edges properly can result in possible bleeding or leakage.
- 3. Inspect the tissue wall to ensure that the forceps do not damage tissue during manipulation.
- 4. Place the instrument jaws onto the everted tissue edges to be anastomosed, making certain that the tissue fits completely within the confines of the jaws. The tissue must comfortably fit within the confines of the jaws, or the use of the instrument is contraindicated.
- 5. Squeeze the levers together fully until a discernible click is felt. As the levers

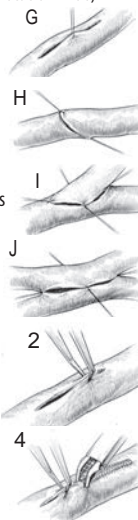are squeezed, the clip is held firmly in the jaws and closed around the tissue. Clip placement should be as close as possible. There should not be more than 0.5 mm between clips (see illustration #7).

**FAILURE TO COMPLETELY SQUEEZE THE LEVERS CAN RESULT IN CLIP MALFORMATION AND POSSIBLE BLEEDING OR LEAKAGE.**

6. Release the levers to disengage the clip from the AnastoClip GC applicator and remove the clip applicator. (The closed clip is disengaged automatically from the jaws.) The applicator automatically advances the next clip for successive applications.
7. Check tightness of clip placement. If desired, the clip can be removed with the AnastoClip Remover and a new clip can be placed with the AnastoClip GC applicator.
8. After completion of the anastomosis, one or more clips, and/or sutures, may be used to control bleeding or leakage from the anastomotic site (if necessary.)
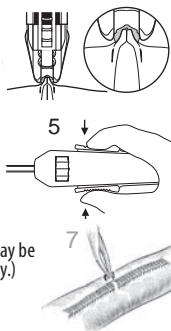
#### A) Size M

Approximate Span Before Closure	Approximate Overall Length	Clips Per Applicator
1.1 mm	2.3 mm	35

The shape of the closed clip may vary according to tissue thickness.

#### B) Size L

Approximate Span Before Closure	Approximate Overall Length	Clips Per Applicator
1.7 mm	3.3 mm	35

The shape of the closed clip may vary according to tissue thickness.

#### C) Size XL

Approximate Span Before Closure	Approximate Overall Length	Clips Per Applicator
2.5 mm	4.93 mm	25

The shape of the closed clip may vary according to tissue thickness.

#### Warning

- SYMMETRICALLY EVERT ALL TISSUE EDGES
- PLACE CLIPS AS CLOSE AS POSSIBLE TO ONE ANOTHER
- AVOID USE ON VESSELS SEVERELY COMPROMISED DURING ENDARTERECTOMY
- FAILURE TO OBSERVE THE ABOVE MAY CONTRIBUTE TO INTEROPERATIVE OR POST-OPERATIVE (SEVERAL DAYS) FAILURE OF ANASTOMOSIS RESULTING IN SERIOUS PATIENT INJURY.

#### Cautions

1. Squeeze the levers together fully until a discernible click is felt. Failure to squeeze the levers completely can result in clip malformation and possible bleeding or leakage.
2. Ensure that the tissue to be anastomosed fits completely within the confines of the jaws or bleeding and leakage may result.
3. Place the clips in such a fashion that they are not "rocking" on their axis ("tips").
4. Inspect the anastomotic site to ensure proper application and that hemostasis has been achieved. If bleeding is observed after application, additional clips or placement of manual sutures may be necessary to complete hemostasis.
5. Do not evert the tissue by grasping one tissue edge with one pair of forceps and the other tissue

edge with another pair of forceps and pulling them together to evert and apply the clips. This may result in asymmetrical eversion of tissue, which could result in possible bleeding or leakage.

6. Inspect the tissue wall to ensure that the forceps do not damage tissue during manipulation.
7. When using the AnastoClip GC applier with tissue or prosthetic graft material (e.g., Meadox, knitted polyester), ensure that the total thickness of the everted tissue to be anastomosed does not exceed the total width of the clip being used (see tables A, B, and C).

### **Contraindications**

1. Do not use the AnastoClip GC applier if tissue can not be properly everted due to the presence of arteriosclerotic or calcified material, or where the vessel has been severely compromised due to endarterectomy (e.g., carotid or any other artery in this condition).
2. The clips should not be more than 0.5 mm between one another. If this cannot be achieved, use of the AnastoClip GC applier is contraindicated.
3. Do not use the AnastoClip GC applier if all the tissue layers cannot be completely symmetrically everted prior to application of the clip.
4. Do not use the AnastoClip GC applier on tissue that is too friable for use of sutures.
5. The tissue must comfortably fit within the confines of the jaws, or use of the instrument is contraindicated.
6. The AnastoClip GC applier is not intended for use except as indicated. **DISCARD AFTER USE. DO NOT RESTERILIZE.**

### **Complications**

- Infection
- Pain/swelling
- Tissue damage
- Blood leakage
- Hemorrhage
- Embolism
- Stenosis
- Intimal dissection
- Thrombosis/occlusion
- Anastomosis rupture
- Dehiscence
- Necrosis
- Seroma
- Anastomotic insufficiency
- Pseudoaneurysm
- Nerve injury
- Steal phenomenon
- Intimal hyperplasia
- Cerebrospinal fluid (CSF) leak
- Pseudomeningocele
- Neuropathic pain
- Spinal infarction

### **Resterilization/Re-use**

This device is single-use only. Do not reuse, reprocess, or re-sterilize. The cleanliness and sterility of the re-processed device cannot be assured. Reuse of the device may lead to cross contamination, infection, or patient death. The performance characteristics of the device may be compromised due to reprocessing or re-sterilization since the device was only designed and tested for single use. The shelf life of the device is based on single use only.

### **Notices: Limited Product Warranty; Limitation Of Remedies**

LeMaitre Vascular, Inc. warrants that reasonable care has been used in the manufacture of this device. Except as explicitly provided herein, LEMAITRE VASCULAR (AS USED IN THIS SECTION, SUCH TERM INCLUDES LEMAITRE VASCULAR, INC., ITS AFFILIATES, AND THEIR RESPECTIVE EMPLOYEES, OFFICERS, DIRECTORS, MANAGERS, AND AGENTS) MAKES NO EXPRESS OR IMPLIED WARRANTIES WITH RESPECT TO THIS DEVICE, WHETHER ARISING BY OPERATION OF LAW OR OTHERWISE (INCLUDING, WITHOUT LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE) AND HEREBY DISCLAIMS THE SAME. LeMaitre Vascular makes no representation regarding the suitability for any particular treatment in which this device is used, which determination is the sole responsibility of the purchaser. This limited warranty does not apply to the extent of any abuse or misuse of, or failure to properly store, this device by the purchaser or any third party. The sole remedy for a breach of this limited warranty shall be replacement of, or refund of the purchase price for, this device (at LeMaitre Vascular's sole option) following the purchaser's return of the device to LeMaitre Vascular. This warranty shall terminate on the expiration date for this device.

IN NO EVENT SHALL LEMAITRE VASCULAR BE LIABLE FOR ANY DIRECT, INDIRECT, CONSEQUENTIAL, SPECIAL, PUNITIVE, OR EXEMPLARY DAMAGES. IN NO EVENT WILL THE AGGREGATE LIABILITY OF LEMAITRE VASCULAR WITH RESPECT TO THIS DEVICE, HOWEVER ARISING, UNDER ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, TORT, STRICT LIABILITY, OR OTHERWISE, EXCEED ONE THOUSAND DOLLARS (US\$1,000), REGARDLESS OF WHETHER LEMAITRE VASCULAR HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH LOSS, AND NOTWITHSTANDING THE FAILURE OF THE ESSENTIAL PURPOSE OF ANY REMEDY. THESE LIMITATIONS APPLY TO ANY THIRD-PARTY CLAIMS.

A revision or issue date for these instructions is included on the back page of these Instructions for Use for the user's information. If twenty-four (24) months has elapsed between this date and product use, the user should contact LeMaitre Vascular to see if additional product information is available.

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LeMaitre, AnastoClip, and AnastoClip GC are registered trademarks of LeMaitre Vascular, Inc. U.S. Patents -- Visit [www.lemaitre.com/patents.asp](http://www.lemaitre.com/patents.asp)
© 2017 LeMaitre Vascular, Inc. All rights reserved. Printed in the USA.
LeMaitre Vascular GmbH