AnastoClip[®] AC Closure System

English — Instructions for Use

AnastoClip[®] AC VerschlussSystem

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STERILE EO Rx only

IMPORTANT!

This booklet is designed to assist in using the AnastoClip® AC Closure System with titanium clips. It is not a reference to surgical stapling techniques.

BEFORE USING PRODUCT, READ THE FOLLOWING INFORMATION THOROUGHLY.

Indications

The AnastoClip applicator has application in the repair of arteriovenous fistulae and the creation of everting anastomoses in blood vessels and other small tubular structures. The Applicator is also intended for approximation of dural tissue (dural approximation is cleared only in the United States, the European Union, Australia and New Zealand).

Effects

The AnastoClip AC applicator is available in three (3) clip sizes: medium-1.1 mm, large-1.7 mm and extra large-2.5 mm. The AnastoClip applicator consists of a rotating shaft and an integral cartridge containing titanium clips. As the levers of the applicator are squeezed together, the clip is closed around the everted tissue edges. As the levers are released, a new clip is automatically loaded into the clip applicator jaws. It is recommended to use the Tissue Everting Forceps to aid in the everting of the tissue edges and the AnastoClip remover for the removal of any clips (if necessary).

MR Compatibility

Non-clinical testing demonstrated that the AnastoClip AC Closure System is MR Conditional. A patient with this device can be scanned safely in an MR system under the following conditions:

Static magnetic field of 1.5-Tesla and 3-Tesla, only


Maximum spatial gradient magnetic field of 2,000-Gauss/cm (20-T/m)

Maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2-W/kg for 15 minutes of scanning (i.e., per pulse sequence) in the Normal Operating Mode

Under the scan conditions defined, the AnastoClip AC Closure System is expected to produce a maximum temperature rise of 1.5°C after 15-minutes of continuous scanning (i.e., per pulse sequence). In non-clinical testing, the image artifact caused by the AnastoClip AC Closure System extends approximately 4-mm from this implant when imaged using a gradient echo pulse sequence and a 3-Tesla MR system.

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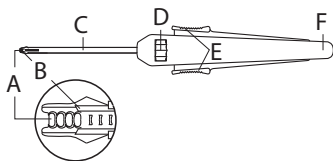How Supplied:

The applicator is supplied sterile. The sterility of the device is assured as long as the packaging is not opened or damaged.

The entire device is considered to be non-pyrogenic for vascular applications.


The implantable clips, the stainless steel shaft and jaws of the applicator are considered to be non-pyrogenic for dura applications. The plastic part of the applicator (handle, levers and rotation knob) should not be in contact with cerebral spinal fluid.

Instructions For Use

NOTE: It is recommended to wear loupes. A 2.5X magnification is suggested.

1. Preparation of tissues is recommended as follows:

- G) ARTERIOTOMY OR VENOTOMY: One optional stay suture at mid-incision.
- H) END-TO-END: Horizontal mattress sutures at 3 and 9 o'clock.
- I) END-TO-SIDE: Horizontal mattress sutures heel and toe: stay sutures at 3 and 9 o'clock.
- J) SIDE-TO-SIDE: Horizontal mattress sutures at 12 and 6 o'clock: stay sutures at 3 and 9 o'clock.

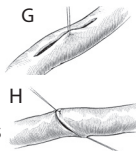NOTE: Additional sutures may be placed depending on the length of the closure, in order to facilitate a symmetrical eversion.

2. Symmetrically evert all tissue layers for secure nonpenetrating clip placement. Evert the tissue edges of the vessel with either design of tissue everting forceps. Ensure that all tissue edges are symmetrically everted prior to applying the clip. Failure to symmetrically evert the tissue edges properly can result in possible bleeding or leakage.

3. Inspect the tissue wall to ensure that the forceps do not damage tissue during manipulation.

NOTE: Atraumatic Tissue Everting Forceps are designed to minimize potential damage to blood vessels or other small tubular structures.

4. Place the instrument jaws onto the everted tissue edges to be anastomosed, making certain that the tissue fits completely within the confines of the jaws. The tissue must comfortably fit within the confines of the jaws, or the use of the instrument is contraindicated.

5. Squeeze the levers together fully until a discernible click is felt. As the levers are squeezed, the clip is held firmly in the jaws and closed around the tissue. Clip placement should be as close as possible. There should not be more than 0.5 mm between clips (see illustration #7).


FAILURE TO COMPLETELY SQUEEZE THE LEVERS CAN RESULT IN CLIP MALFORMATION AND POSSIBLE BLEEDING OR LEAKAGE.

6. Release the levers to disengage the clip from the AnastoClip applicator and remove the clip applicator. (The closed clip is disengaged automatically from the jaws.) The applicator automatically advances the next clip for successive applications.

7. Check tightness of clip placement. Tissue should completely fill clip opening and clip should not loosely rock side to side.

8. If desired, the clip can be removed with the AnastoClip Remover and a new clip can be placed with the AnastoClip AC applicator.

9. After completion of the anastomosis, one or more clips, and/or sutures, may be used to control bleeding or leakage from the anastomotic site (if necessary).

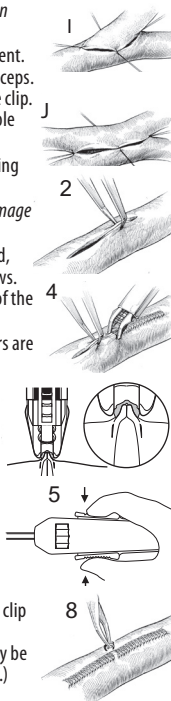A) 4012-01: Size M


Approximate Span Before Closure	Approximate Overall Length	Clips Per Applicator
1.4 mm	2.3 mm	35

The shape of the closed clip may vary according to tissue thickness.

The amount of Titanium that a patient is exposed to is 46 milligrams per Applicator if the entire load of the applicator is implanted.

B) 4012-02: Size L

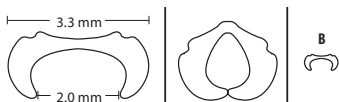
Approximate Span Before Closure	Approximate Overall Length	Clips Per Applicator
2.0 mm	3.3 mm	35

The shape of the closed clip may vary according to tissue thickness.

The amount of Titanium that a patient is exposed to is 97 milligrams per Applicator if the entire load of the applicator is implanted.

C) 4012-03: Size XL

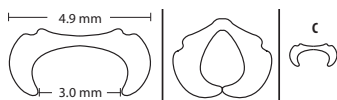Approximate Span Before Closure	Approximate Overall Length	Clips Per Applier
3.0 mm	4.9 mm	25

The shape of the closed clip may vary according to tissue thickness.

The amount of Titanium that a patient is exposed to is 140 milligrams per Applier if the entire load of the applier is implanted.

WARNING

- **SYMMETRICALLY EVERT ALL TISSUE EDGES**
 - **PLACE CLIPS AS CLOSE AS POSSIBLE TO ONE ANOTHER**
 - **AVOID USE ON VESSELS SEVERELY COMPROMISED DURING ENDARTERECTOMY**
- FAILURE TO OBSERVE THE ABOVE MAY CONTRIBUTE TO INTEROPERATIVE OR POST-OPERATIVE (SEVERAL DAYS) FAILURE OF ANASTOMOSIS RESULTING IN SERIOUS PATIENT INJURY.**

Cautions

1. Squeeze the levers together fully until a discernible click is felt. Failure to squeeze the levers completely can result in clip malformation and possible bleeding or leakage.
2. Ensure that the tissue to be anastomosed fits completely within the confines of the jaws or bleeding and leakage may result. See figure below.
3. Place the AnastoClips in such a fashion that they are not “rocking” on their axis (“tips”). See figure below.

4. Inspect the anastomotic site to ensure proper application and that hemostasis has been achieved. If bleeding is observed after application, additional clips or placement of manual sutures may be necessary to complete hemostasis.
5. Do not evert the tissue by grasping one tissue edge with one pair of forceps and the other tissue edge with another pair of forceps and pulling them together to evert and apply the clips. This may result in asymmetrical eversion of tissue, which could result in possible bleeding or leakage.
6. Inspect the tissue wall to ensure that the forceps do not damage tissue during manipulation.
7. When using the AnastoClip applier with tissue or prosthetic graft material (e.g., Meadox, knitted polyester) ensure that the total thickness of the everted tissue and graft material to be anastomosed does not exceed 1.0 mm for the medium 1.4 mm clip applier, 1.6 mm for the large 2.0 mm clip applier and 1.8 mm for the extra large 3.0 mm clip applier.

Contraindications

1. Do not use the AnastoClip applier if tissue can not be properly everted due to the presence of arteriosclerotic or calcified material, or where the vessel has been severely compromised due to endarterectomy (e.g., carotid or any other artery in this condition).
2. The AnastoClips should not be more than 0.5 mm between one another. If this cannot be achieved, use of the AnastoClip applier is contraindicated.
3. Do not use the AnastoClip applier if all the tissue layers cannot be completely symmetrically everted prior to application of the clip.
4. Do not use the AnastoClip applier on tissue that is too friable for use of sutures.
5. The tissue must comfortably fit within the confines of the jaws, or use of the instrument is contraindicated.

Complications

- Infection
- Pain/swelling
- Tissue damage
- Blood leakage
- Hemorrhage
- Embolism
- Stenosis
- Intimal dissection
- Thrombosis/occlusion
- Anastomosis rupture
- Dehiscence
- Necrosis
- Seroma
- Anastomotic insufficiency
- Pseudoaneurysm
- Nerve injury
- Steal phenomenon
- Intimal hyperplasia
- Cerebrospinal fluid (CSF) leak
- Pseudomeningocele

Resterilization/Re-use

This device is single-use only. Do not reuse, reprocess, or re-sterilize. The cleanliness and sterility of the re-processed device cannot be assured. Reuse of the device may lead to cross contamination, infection, or patient death. The performance characteristics of the device may be compromised due to reprocessing or re-sterilization since the device was only designed and tested for single use. The shelf life of the device is based on single use only.

Safe Handling and Disposal

If serious medical incidents should arise during use of this medical device, users should notify both LeMaitre Vascular and the Competent Authority of the country where the user is located. This product contains no sharps, heavy metals or radioisotopes, and is not infectious or pathogenic. No special requirements for disposal are evident. Please consult local regulations to verify proper disposal.

Packaging and Shipping of Explanted AnastoClips:

Return of the shipment to LeMaitre Vascular depends on 2 crucial questions:

1. Is the explant obtained from a patient with a known or presumed pathogenic condition at the time of the explant?
2. Is the explant obtained from a patient that has a known treatment history that involves therapeutic radionuclides within the last 6 months?

In the event that the response to question 1 or 2 is affirmative, LeMaitre Vascular does not provide adequate guidance for shipment. THESE EXPLANTS SHOULD NOT BE RETURNED TO LEMAITRE VASCULAR UNDER ANY CIRCUMSTANCES. In these cases, the explant should be disposed according to local regulations.

For explants that do not have pathogenic or radiologic hazards, please use the following:

Pre-explantation:

LeMaitre Vascular can accept clinical information that is patient-anonymized.

LeMaitre Vascular requests information including:

1. The original diagnosis which resulted in use of the implant,
2. The patient's medical history relevant to the implant, including the hospital or clinic at which the device was implanted.
3. The patient's implant experience prior to implant removal.
4. The hospital or clinic at which the explantation was done, and date of retrieval.

Explantation:

1. Explanted AnastoClips should first be cleaned using one of the following:
 - a. Intense water rinse, followed by transfer to 70-80% ethanol, or
 - b. Cleaning in isopropanol with subsequent ultrasonic treatment, or
 - c. Cleaning in proteolytic enzyme solution, or
 - d. Cleaning in sodium hypochlorite solution (50-60 mg/l), or
 - e. Cleaning in 3% hydrogen peroxide.
2. Explanted AnastoClips should then be decontaminated with steam autoclaving or ethylene oxide.
3. AnastoClips should be packaged dry.

Packaging:

1. Explants should be sealed and packed in a manner that minimizes potential for breakage, contamination of the environment or exposure to those handling such packages during transit. The sealed primary container should be placed inside secondary packaging. Primary and secondary packaging must then be packaged inside an outer package, which must be a rigid, fiberboard box.
2. The secondary packaging of cleaned, decontaminated AnastoClip explants should be labeled with an ISO 7000-0659 Biohazard symbol. Outer packaging should also be labelled with "Exempt human specimen" printed on a label with dimensions at least 10 x 10 cm.
3. Shipping paper for the outer shipping container is not required.
4. Packages prepared in the above manner may be shipped to:

LeMaitre Vascular
63 Second Avenue
Burlington, MA 01803

Notices: Limited Product Warranty; Limitation Of Remedies

LeMaitre Vascular, Inc. warrants that reasonable care has been used in the manufacture of this device and that this device is suitable for the indication(s) expressly specified in these instructions for use. Except as explicitly provided herein, LEMAITRE VASCULAR (AS USED IN THIS SECTION, SUCH TERM INCLUDES LEMAITRE VASCULAR, INC., ITS AFFILIATES, AND THEIR RESPECTIVE EMPLOYEES, OFFICERS, DIRECTORS, MANAGERS, AND AGENTS) MAKES NO EXPRESS OR IMPLIED WARRANTIES WITH RESPECT TO THIS DEVICE, WHETHER ARISING BY OPERATION OF LAW OR OTHERWISE (INCLUDING, WITHOUT

LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE) AND HEREBY DISCLAIMS THE SAME. This limited warranty does not apply to the extent of any abuse or misuse of, or failure to properly store, this device by the purchaser or any third party. The sole remedy for a breach of this limited warranty shall be replacement of, or refund of the purchase price for, this device (at LeMaitre Vascular's sole option) following the purchaser's return of the device to LeMaitre Vascular. This warranty shall terminate on the expiration date for this device.

IN NO EVENT SHALL LEMAITRE VASCULAR BE LIABLE FOR ANY DIRECT, INDIRECT, CONSEQUENTIAL, SPECIAL, PUNITIVE, OR EXEMPLARY DAMAGES. IN NO EVENT WILL THE AGGREGATE LIABILITY OF LEMAITRE VASCULAR WITH RESPECT TO THIS DEVICE, HOWEVER ARISING, UNDER ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, TORT, STRICT LIABILITY, OR OTHERWISE, EXCEED ONE THOUSAND DOLLARS (US\$1,000), REGARDLESS OF WHETHER LEMAITRE VASCULAR HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH LOSS, AND NOTWITHSTANDING THE FAILURE OF THE ESSENTIAL PURPOSE OF ANY REMEDY. THESE LIMITATIONS APPLY TO ANY THIRD-PARTY CLAIMS.

These limitations do not apply to consumers in Australia or to the extent they are precluded by local law in any other jurisdiction.

A revision or issue date for these instructions is included on the back page of these Instructions for Use for the user's information. If twenty-four (24) months has elapsed between this date and product use, the user should contact LeMaitre Vascular to see if additional product information is available.

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LeMaitre and AnastoClip are registered trademarks of LeMaitre Vascular, Inc.
© 2020 LeMaitre Vascular, Inc. All rights reserved. Printed in the USA.
LeMaitre Vascular GmbH