**Over-the-Wire LeMaitre® Valvulotome**

Instructions for Use - English

**Over-the-Wire LeMaitre® Valvulotom**

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**Valvulotome Over-the-Wire LeMaitre®**

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**Valvulotomo over-the-wire (OTW) LeMaitre®**

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**Valvulotomo guiado LeMaitre®**

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**Valvulótomo através de um fio guia LeMaitre®**

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**STERILE**

### **Description**

The Over-the-Wire LeMaitre® Valvulotome is used to cut saphenous vein valves during bypass procedures. The Centering Hoops keep the head of the Over-the-Wire LeMaitre Valvulotome centered in the vessel and prevents the valve-cutting blades from damaging the vessel wall. The size of the Centering Hoops and Cutting Blades adjust to the internal diameter of the vein as the Over-the-Wire LeMaitre Valvulotome is being drawn through the vessel. The LeMills Valvulotome is provided as an accessory device to cut the last valve set as needed. Kit includes LeMills Valvulotome.

### **Indication for Use**

The Over-the-Wire LeMaitre Valvulotome is indicated for cutting saphenous vein valves during in-situ bypass. The LeMills Valvulotome is intended to cut venous valves.

### **Contraindications**

1. pregnant women,
2. patients with allergy to any of the device materials,
3. endarterectomy procedures,
4. thrombolysis procedures,
5. vein stripping procedures,
6. embolectomy procedures, or
7. vessel dilation procedures.

### **Warnings**

1. Do not use if the package or device is damaged.
2. Do not use if packaging is opened outside a sterile environment.
3. The Over-the-Wire Valvulotome and the LeMills Valvulotome are to be used by a qualified physician.
4. Do not insert the Over-the-Wire Valvulotome into a vessel or extract from a vessel in the open position.
5. The Over-the-Wire Valvulotome must be used with an 0.035" Guidewire.
6. Do not open or close the Over-the-Wire Valvulotome while in a coiled configuration.
7. The Over-the-Wire Valvulotome must be flushed with saline or heparinized saline only.
8. Do not pass the Over-the-Wire Valvulotome or LeMills Valvulotome through a vessel that has undergone synthetic grafting or contains implants.
9. Do not rotate the Over-the-Wire Valvulotome in a vessel.
10. A maximum of 8 passes with the Over-the-Wire Valvulotome may be completed.

### **Precautions**

1. United States Federal and other law restricts this device to sale on or by the order of a physician.
2. Do not use device past the expiration date printed on the labeling.
3. These devices are single use only. Do not reuse, reprocess, or re-sterilize (reference "Re-sterilization/ Repackaging" section of this document for further details).
4. Keep any dissected portion of the vein straight and avoid twisting the vein during valvulotomy.
5. Use caution when advancing the guidewire past the proximal anastomosis.
6. Use caution when tracking the Over-the-Wire Valvulotome in thrombophlebitic veins.
7. Do not advance the Over-the-Wire Valvulotome with the blades in the open position.

### **Potential Complications**

- Vessel damage
- Vessel occlusion/ stenosis
- Hematoma


- Hemorrhage
- Thrombus formation
- Wound infections
- Erythema
- Myocardial infarction
- Amputation
- Death

**Additional Required Items**


1. .035" Guidewire
2. Radiopaque Tape

**Specifications**

Over-the-Wire Valvulotome:

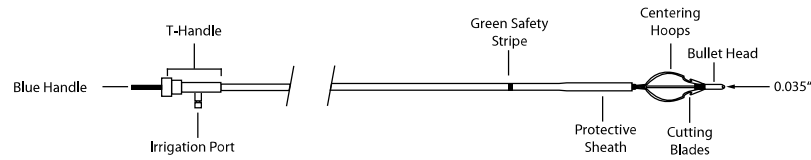LeMills Valvulotome:

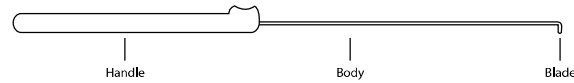Catalog Number	Useable Length	Closed Catheter Outer Diameter	Maximum Open Cutting Blade Diameter	Maximum Open Centering Hoop Diameter
1005-00, e1005-00	100 cm	2.0 mm	6.0 mm	9.5 mm
1006-00, e1006-00	40 cm	2.0 mm	6.0 mm	9.5 mm

**To Open Package**

1. Open box and remove sealed sterile unit.
2. Open sealed sterile unit and present contents to personnel in sterile environment.
3. Carefully remove lids, as applicable.

**Caution:** Remove product from tray with care to ensure centering hoops are not damaged.

**Pre-Use Check**

Inspect centering hoops and cutting blades for damage and alignment prior to use.

**Instructions for Use**

1. Perform an intra-operative venogram using Radiopaque Tape; mark the fistulae.
2. Excise the Sapheno-Femoral Complex.
3. Excise first valve set under direct vision using standard technique.
4. Perform the proximal anastomosis.
5. Ensure that the leg is fully extended and pass the 0.035" guidewire well beyond the treatment area.
6. Flush the Over-the-Wire Valvulotome (in the open position) by removing Irrigation Port cap, blocking guidewire lumen on the Blue Handle, and flushing with sterile heparinized saline. Replace cap upon completion.

7. Sheath the Centering Hoops within the Protective Sheath by advancing the T-Handle until Cutting Blades are fully enclosed and a positive stop is detected.
8. Load the device over the guidewire carefully.
9. Advance into the vein to the treatment area; 2 cm to 3 cm distal to the anastomosis.
10. Expose the Cutting Blades and Centering Hoops by maintaining the position of the Blue Handle and retracting the T-Handle until a positive stop is detected.  
**Note:** *Failure to maintain the position of the Blue Handle while exposing cutting blades may cause damage to the anastomosis.*
11. Slowly retract the Over-the-Wire Valvulotome in the open position to disrupt valve sets.  
**Warning:** *Do not rotate the Over-the-Wire Valvulotome.*
12. When the Green Safety Stripe is visible, sheath the Cutting Blades and Centering Hoops into the Protective Sheath by maintaining the position of the Blue Handle and advancing the T-Handle.  
**Warning:** *Failure to sheath Cutting Blades and Centering Hoops prior to retracting the Over-the-Wire Valvulotome from the vessel may cause damage.*  
**Note:** *If the Over-the-Wire Valvulotome is set aside for further use within the same procedure, device must be stored in the open position in a basin of heparinized saline. Flush device as noted in step 6 prior to reintroduction.*
13. Ligate the venous tributaries previously located and marked.  
**Note:** *Do not ligate the last major tributary until after the distal anastomosis is complete in order to allow blood run off.*
14. Pass the device through the vessel an additional 1-2 times after ligating the venous tributaries.
15. Steps 9 through 12 may be repeated as necessary. Use the LeMills Valvulotome to excise last valve set following standard technique as needed.  
**Warning:** *A maximum of 8 passes with the Over-the-Wire Valvulotome may be completed.*
16. Confirm free flow through saphenous vein using doppler.
17. Complete distal anastomosis to the appropriate arterial segment.
18. Perform a completion arteriogram using radiopaque tape to locate defects.

#### **Storage**

Store in a cool, dry place.

#### **Re-sterilization/Repackaging**

These devices are single-use only. Do not reuse, reprocess, or re-sterilize. The cleanliness and sterility of the re-processed device cannot be assured. Reuse of the device may lead to cross contamination, infection, or patient death. The performance characteristics of the device may be compromised due to reprocessing or re-sterilization since the device was only designed and tested for single use. The shelf life of the device is based on single use only. If for any reason this device must be returned to LeMaitre Vascular, place it in its original packaging and return it to the address listed on the box.

#### **Limited Product Warranty; Limitation of Remedies**

LeMaitre Vascular, Inc. warrants that reasonable care has been used in the manufacture of this device. Except as explicitly provided herein, LEMAITRE VASCULAR (AS USED IN THIS SECTION, SUCH TERMS INCLUDE LEMAITRE VASCULAR, INC., ITS AFFILIATES, AND THEIR RESPECTIVE EMPLOYEES, OFFICERS, DIRECTORS, MANAGERS, AND AGENTS) MAKES NO EXPRESS OR IMPLIED WARRANTIES WITH RESPECT TO THIS DEVICE, WHETHER ARISING BY OPERATION OF LAW OR OTHERWISE (INCLUDING, WITHOUT LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE) AND HEREBY DISCLAIMS THE SAME. LeMaitre Vascular makes no representation regarding the suitability for any particular treatment in which this device is used, which determination is the sole responsibility of the purchaser. This limited warranty does not apply to the extent of any abuse or misuse of, or failure to properly store, this device by the purchaser or any third party. The sole remedy for a breach of this limited warranty shall be replacement of, or refund of the purchase price for, this device (at LeMaitre Vascular's sole option) following the purchaser's return of the device to LeMaitre Vascular. This warranty shall terminate on the expiration date for this device.

IN NO EVENT SHALL LEMAITRE VASCULAR BE LIABLE FOR ANY DIRECT, INDIRECT, CONSEQUENTIAL, SPECIAL, PUNITIVE, OR EXEMPLARY DAMAGES. IN NO EVENT WILL THE AGGREGATE LIABILITY OF LEMAITRE VASCULAR WITH RESPECT TO THIS DEVICE, HOWEVER ARISING, UNDER ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, TORT, STRICT LIABILITY, OR OTHERWISE, EXCEED ONE THOUSAND DOLLARS (US\$1,000), REGARDLESS OF WHETHER LEMAITRE VASCULAR HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH LOSS, AND NOTWITHSTANDING THE FAILURE OF THE ESSENTIAL PURPOSE OF ANY REMEDY. THESE LIMITATIONS APPLY TO ANY THIRD-PARTY CLAIMS.

A revision or issue date for these instructions is included on the back page of these Instructions for Use for the user's information. If twenty-four (24) months has elapsed between this date and product use, the user should contact LeMaitre Vascular to see if additional product information is available.

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LeMaitre is a registered trademark of LeMaitre Vascular, Inc.
©2013 LeMaitre Vascular, Inc. All rights reserved. Printed in the USA.
LeMaitre Vascular GmbH